



Application Number


011001

200 Series GC
Application Note



Red Wine

Column	Alltech EC-Wax 30m x 0.45id 1ul film thickness	1 Acetaldehyde
Flow	Hydrogen 4psi	2 Acet/EtAc
Program	50C for 3min 10C/min to 220C Hold 5min	3 Methanol
Injector Temp	220C Split mode	4 Ethanol
FID Detector Temp	240C	5 N-Propanol
		6 iso-Butanol
		7 n-Butanol
		8 iso-Amyl

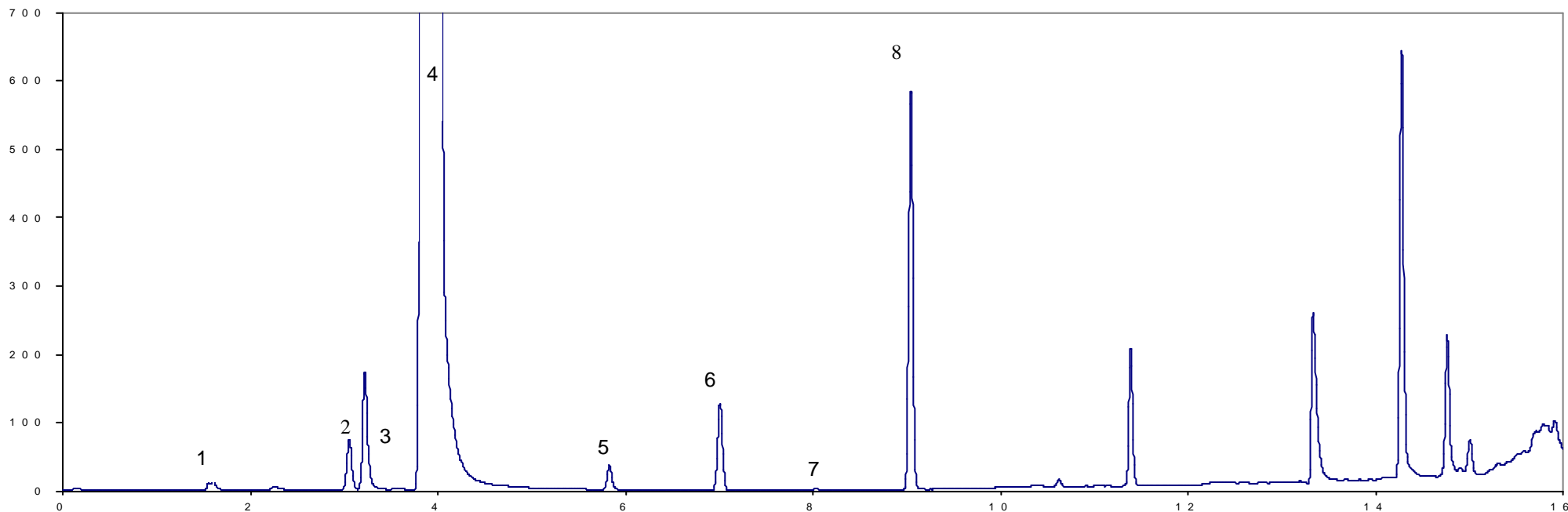

Application Number

011001